

Ferret CD56 Antibody, Biotin

Mouse monoclonal

Clone: FT56-4B3

Catalog: M2039b

Pack Size: 100µg

CD56, [UniProt # A0A8U0RZ03](#), also known as NCAM, is a glycoprotein of the Ig superfamily. It serves as a marker for natural killer cells.

Description	
Immunogen	HEK-293-derived recombinant CD86. Leu20-Ile713
Reactivity	Ferret
Source	Mouse monoclonal IgG ₁ coupled to biotin
Purification	Protein G purified from hybridoma cell culture supernatant
Applications	Flow cytometry
Formulation	Lyophilized from PBS, 0.1% BSA, 0.05% ProClin300 with Trehalose Reconstitute in 100µL ddH ₂ O to 1 mg/mL
Shipping	Ambient
Storage	5 years at -20°C to -80°C as supplied 1 month at 4°C after reconstitution with preservative 1 year at -20°C to -80°C after reconstitution
Expiration	2 years from date of receipt


For Research Use Only

AMPERSAND BIOSCIENCES, LLC

6459 State Route 30, Lake Clear, NY, 12945 | PHONE: (518) 891 0204 |

Rev.1

Flow cytometric analysis of CD56 expression on ferret PBMCs.

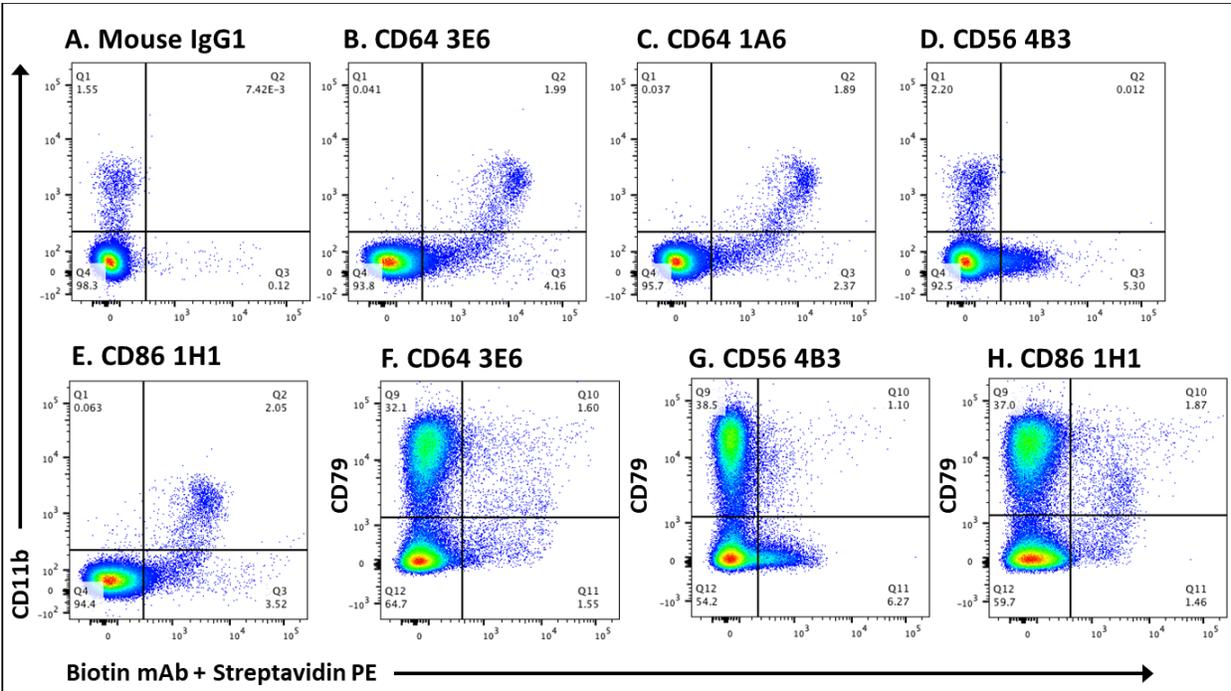

Representative double staining for female ferret PBMCs stained with biotinylated mouse anti-ferret antibodies followed by streptavidin PE) in conjunction with ferret cross-reactive antibodies described in the literature. Cells were first gated on a live/dead stain and doublets removed before gating on forward and side scatter. Shown is the total lymphocyte population double staining with anti-mouse CD11b and each mouse biotinylated antibody (Fig A-E). CD64 (*Cat No. M2038*) and CD86 (*Cat No. M2040*) co-stain with macrophage marker CD11b as expected, while the NK marker CD56 does not co-stain with CD11b. In panels F-H, cells were double stained with human anti-CD79, a B cell marker. CD64 and CD86 do co-stain to some extent with the B cell marker as B cells are known to express FcγRI and the co-stimulatory molecule CD86. A small percentage of CD56+ cells also stain positive for CD79 marker.

For Research Use Only

Ampersand Biosciences products may not be sequenced, reverse engineered, or investigated in any way to determine composition without the expressed written consent of Ampersand Biosciences.