

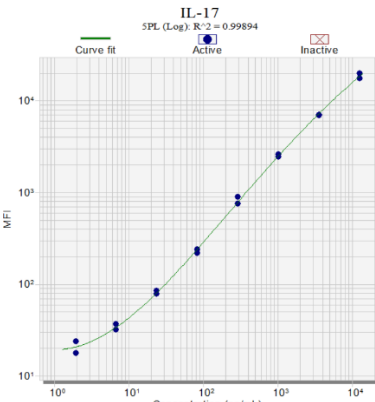
Hamster IL-17 Antibody

Mouse monoclonal, clone: HM17-3E4

Catalog: M2064

Pack Size: 100µg

Interleukin 17 (IL-17) is a member of cytokines that promote immunity against infections from pathogens. IL-17 has also been shown to drive inflammation during autoimmune responses. This makes it a target for drugs to treat autoimmune and other inflammatory conditions.

Description	
Immunogen	HEK293 derived recombinant hamster IL-17. Met1-Ala158 with C-terminal His tag
Reactivity	Hamster
Source	Mouse monoclonal IgG ₁
Purification	Protein G purified from hybridoma cell culture supernatant
Applications	ELISA (Detection)
Formulation	Lyophilized from PBS with Trehalose Reconstitute in 100µL ddH ₂ O to 1 mg/mL
Shipping	Ambient
Storage	2 years at -20°C to -80°C as supplied 1 month at 4°C after reconstitution with preservative 1 year at -20°C to -80°C after reconstitution
Expiration	2 years from date of receipt
 IL-17 SPL (Log): R² = 0.99894 Curve fit Active Inactive MFI Concentration (pg/mL)	
Hamster IL-17 Luminex Assay Capture: Hamster IL-17 clone HM17-2H5-1A3 (Cat # M2063) Detection: Hamster IL-17 clone HM17-3E4 (Cat #M2064) Antigen: Recombinant Hamster IL-17 (Cat # P4051)	

For Research Use Only

Ampersand Biosciences products may not be sequenced, reverse engineered or investigated in any way to determine composition without the expressed written consent of Ampersand Biosciences.

AMPERSAND BIOSCIENCES LLC

6459 State Route 30. Lake Clear, NY 12945 | PHONE: (518) 891 0204 |

Rev. 1.0