



Hamster CD69 Antibody, Biotin

Mouse monoclonal, clone: 2B3


Catalog: M2074

Pack Size: 100µg

Cluster of Differentiation 69 (CD69) a transmembrane C-type lectin protein is an early activation marker in the immune system, playing a crucial role in lymphocyte activation, proliferation and retention in lymphoid organs. It is primarily expressed on the surface of activated lymphocytes, including T cells, B cells and NK cells.

Description	
Immunogen	Recombinant hamster CD69 Ala83-Arg206 w/ C-terminal His tag expressed in HEK293
Reactivity	Hamster
Source	Mouse monoclonal IgG ₁ coupled to biotin
Purification	Protein G purified from hybridoma cell culture supernatant
Applications	Flow Cytometry
Formulation	Lyophilized from PBS with Trehalose Reconstitute in 100µL ddH ₂ O to 1 mg/mL
Shipping	Ambient
Storage	2 years at -20°C to -80°C as supplied 1 month at 4°C after reconstitution with preservative 1 year at -20°C to -80°C after reconstitution
Expiration	2 years from date of receipt

Flow cytometric analysis of CD25, CD69 expression on hamster spleen cells

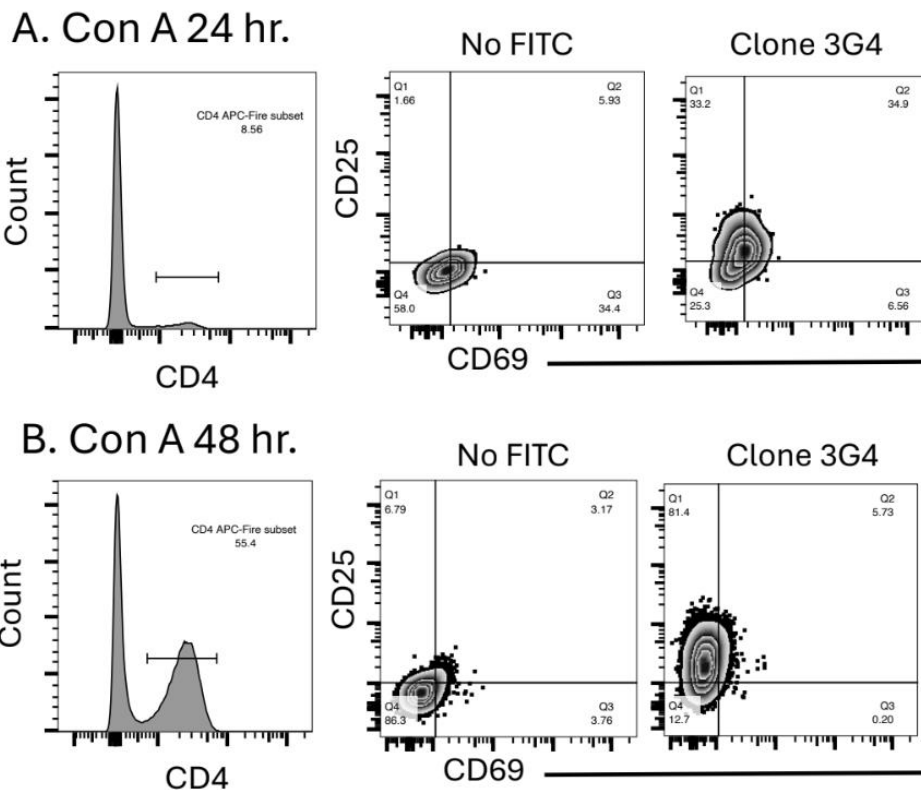

Anti-CD25 clone HM25-3G4 binding increased on CD4 cells after 24 h and 48 h activation with Con A. Anti-CD69 clone HM69-2B3 showed increased levels after 24 h, but less expression after 48 h which was expected with CD69 being an early activation marker. Hamster spleen cells were incubated with Concanavalin A for 24 (A) or 48 (B) hours. For all samples, dead cells were first gated out based on a live-dead dye and live cells further gated on CD4+ cells. CD4+ cells were then analyzed for CD69 and CD25 expression.

For Research Use Only

Ampersand Biosciences products may not be sequenced, reverse engineered, or investigated in any way to determine composition without the expressed written consent of Ampersand Biosciences.

AMPERSAND BIOSCIENCES, LLC

3 Main St., Saranac Lake, NY 12983 | PHONE: (518) 891 0204 |